



INTELLI-MIXER RM-2

INTELLI-MIXER RM-2 is a combination of variable speed rotator and advanced multiple tube mixer and vortex in one machine. Vast selection of racks and operation regimes allow to accommodate the mixer for virtually any kind of mixing from gentle low speed hybridization up to strong test tube vortexing in a whole range of tubes from 0.5ml up to 50ml. Superior performance and flexibility allows to mix almost any kind of laboratory samples and liquids with amazing efficiency. It also allows to program up to 3 personal custom regimes by user to insure always the best up to date functionality.

PRODUCT FEATURES

- ▶ 3 Models selection.
- ▶ Easy to operate interface.
- ▶ Supports whole range of test tubes from 0.5ml up to 50ml.
- ▶ Rotation speed range from 1 - 99 RPM.
- ▶ 18 different programs of rotation and rocking.
- ▶ 2 variable intensity vortexing programs.
- ▶ 3 additional custom programs that could be created by user.
- ▶ FLASH memory on-board for storing personal settings and saving new programs.
- ▶ Maintenance-free high performance stepper motor.
- ▶ Standby low energy consumption mode SLEEP.



PRODUCT BENEFITS

INTELLI-MIX is easy to accommodate the current lab requirements by attaching different racks and selecting appropriate operation regimes. All the adjustments could be made very easily and within seconds. It makes INTELLI-MIXER a perfect assistant in medicine, biology, chemistry and vast variety of different laboratories.

MODELS



RM-2S
L=315mm



RM-2M
L=340mm



RM-2L
L=420mm

ACCESSORIES



Up to 12x 30mm places
Bloating and
Hybridization rack.



Up to 14x 30mm
places



Up to 22x 16mm
places



Up to 38x 13mm
places



Up to 42x 11mm
places



11, 13, 16, 30 (mm)
combined rack



INTELLI-MIXER RM-2

SPECIFICATIONS

Rotation speed	RPM	from 1 up to 99
Rotation speed settings step	RPM	1
Overall number of operation regimes		23
Environmental temperature	°C	from 2 up to 50
Relative air humidity (at 20 °C)	%	up to 80
Power supply adaptor		DC in 12V, 1000 mA


ANY KIND OF TEST TUBES
from 0.5 ml up to 50 ml



18 GREAT MIXING
REGIMES SELECTION.



2 REGIMES OF
ULTRA EFFECTIVE
VORTEX

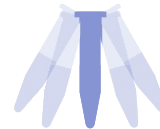

3 BONUS CUSTOM
PROGRAMMABLE
REGIMES

FLASH
MEMORY

Rose Scientific Ltd.
enhancing the sciences...

4027 - 97 Street
Edmonton, Alberta, Canada. T6E 5Y5
Tel: 780-438-5110, 800-661-9289
Fax: 780-462-5776, 800-570-6067

Mississauga Office:
6581 Kitimat Road, Unit #7
Mississauga, Ontario, Canada. L5N 3T5
Tel: 877-826-6706
Fax: 905-814-1978
sales@rosesci.com

www.rosesci.com